MANUAL DE INSTRUCCIONES

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VIBFIT

Manual de instrucciones

ES

Instruction manual

EN

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FR

Manuale D'Istruzioni

IT

Anweisungen

DE

Instruções

PT

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NL

Bruksanvisning

SWE

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PL

Brugsanvisning

DK

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TO KEEP FOR FUTURE REFERENCE.

1. WARNINGS

1. Keep the instructions for future reference. If the device is handed over to a third party, these instructions must also be handed over.
2. Do not use directly on the skin as it may cause friction rubbing.
3. Discontinue use immediately if you feel discomfort.
4. Do not use the product before going to sleep. It may have a stimulating effect and delay sleep.
5. To avoid short-circuit, connect the device correctly and keep it away from metal, rubbish and water.
6. Use the device on a flat surface. If it is not stable, it may tip over and cause injury.
7. Unplug the charger from the wall socket when not in use.
8. If you suspect a problem, unplug the device immediately.
9. Do not place objects on top of the Vibfit.
10. If the connection cable is damaged, it must be replaced. This should only be done by an authorised and qualified person.
11. To avoid the risk of electric shock, do not use the device with wet skin.
12. Keep the product away from moisture and sun exposure. Do not use if wet or damp parts are present.
13. Keep the device out of the reach of children and animals.
14. To avoid injury or irreparable damage, always remain seated when using the Vibfit. Never stand on the unit.
15. Do not attempt to disassemble or repair the appliance. This should only be

done by qualified technicians.

16. For delicate floor surfaces (e.g. parquet), we recommend placing a mat or matting under the unit to prevent damage or discolouration.
17. For indoor use only.
18. The device is designed to be used on the soles of the feet, legs and back. Do not use on the waist, neck or other parts of the body not described in this manual.
19. Do not use this device for any other purpose not specified in this manual.


2. COMMISSIONING AND USE

1. Connect the charger to the socket and the cable to the unit. As soon as it is plugged in, you will hear a beep and the display will light up.
2. Once plugged in, turn it so that the control panel is facing forward. Select a mode according to your needs using the remote control or the buttons on the appliance. Your settings will be shown on the display. To start, press the Start/Stop button on the remote control of the device.
3. To properly turn off the Vibfit, you must first have the device by pressing the Start/Stop button on the device or on the remote control. Then press the power button on the device to turn it off completely.
4. Unplug the appliance when you have finished using it.

3. CORRECT POSTURE FOR USE

Place the sole of your foot or calf on the device. For use of the device in the back area, place it at the level of your back and lean on it.

Note: Never put your full weight on the device, as this may damage it.

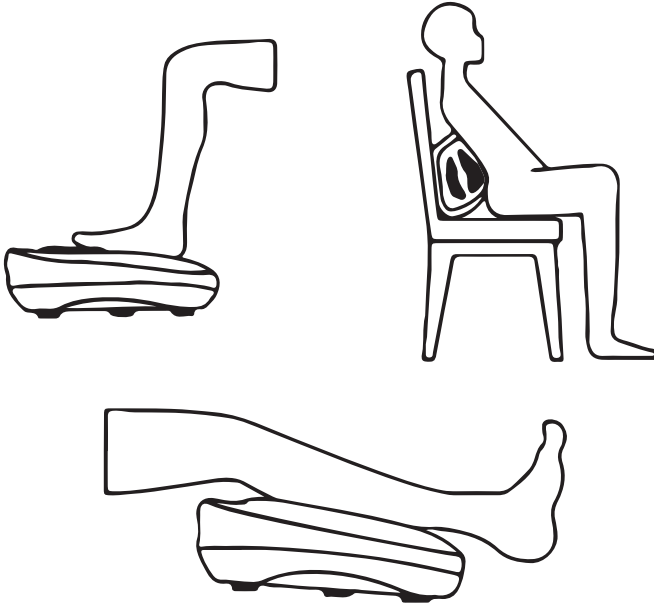4. CONTROL PANEL AND REMOTE CONTROL

- **Power button (Power):** Press this button on your unit to turn it on and off.
- **Start/Stop:** Press this button to start with your selected settings. When the device is running, press this button to stop it.
- **Time - :** Lower the time. The default time is 15 minutes, press this button to reduce it.
- **Time +:** Increase time. Press this button to increase the time. The maximum time setting is 30 minutes.
- **Speed +:** Increases the speed. Speed can be adjusted.
- **Speed - :** Decreases the speed. The speed can be manually adjusted from 1 to 10. Press the button to decrease the speed while in operation.

- **Programme button (P):** You can switch the programmes from manual mode to automatic mode with this button. If you press it repeatedly, you can switch between HO, UI, P1, P2 and P3 programmes. By pressing the programme button once, "UI" will appear on the display. In this mode, you will have the oscillation function only and customisable speed. In programme "P1", the oscillation and massage function are activated and the device executes a programme with a preset range of speeds. It is not possible to change the speed with this programme. In programme "P2", the device will alternate between oscillation and massage with the preset time and speeds. In programme "P3", the acupressure heads are deactivated and only the oscillation remains active, with the preset speed. Manual mode "HO", in standby mode, when you press the "Start/Stop" button, the device will switch to manual mode, "HO" will appear on the display. In this mode, the oscillation and massage function works continuously, the speed is customisable.
- In standby mode, you can adjust the time but not the speed. Once the device has started, you can adjust the speed but not the time.
- By pressing the "P" button on the device several times, you can also switch to manual "HO" mode.
- You can select the "P" (Auto Program) mode after switching on the device, press "Start/Stop" to activate it.
- Stop the Vibfit by first pressing the "Start/Stop" button. Then turn the unit off with the power button or press "Start/Stop" to continue the massage.
- In the "P3" programme, the rotating acupressure heads will stop as soon as you place your feet on the device. This is not a defect.

5. PROGRAMMES P1, P2, P3 AND MANUAL MODE

AUTOMATIC PROGRAMME P1	
Function	Speed range (Changes per minute)
Oscillation and acupressure	1-2-3-4-5-6-7-8-5-3-1-5-8-3

AUTOMATIC PROGRAMME P2		
Standard time	Function	Speed level
Minute 1	Oscillation and acupressure	Level 1
Minute 2	Oscillation only	Level 3
Minute 3	Oscillation and acupressure	Level 5
Minute 4	Oscillation only	Level 7
Minute 5	Oscillation and acupressure	Level 9
Minute 6	Oscillation only	Level 10
Minute 7	Oscillation and acupressure	Level 8

Minute 8	Oscillation only	Level 6
Minute 9	Oscillation and acupressure	Level 4
Minute 10	Oscillation only	Level 2
Minute 11	Oscillation and acupressure	Level 1
Minute 12	Oscillation only	Level 4
Minute 13	Oscillation and acupressure	Level 8
Minute 14	Oscillation only	Level 6
Minute 15	Oscillation and acupressure	Level 4

AUTOMATIC PROGRAMME P3

Function	Speed range (Changes per minute)
Oscillation only	3-5-7-9-2-4-6-8-10-1-5-10-5-1-5

MANUAL MODE "HO"

Function	Speed range
Oscillation and acupressure	Speed is customisable

MANUAL MODE "UI"

Function	Speed range
Oscillation only	Speed is customisable

6. MAINTENANCE

- Always disconnect the plug before cleaning.
- Clean the device and the silicone bases with a damp cloth and neutral soap. Do not use abrasive or corrosive cleaning agents. The silicone pads can be removed for easy cleaning.
- Store the Vibfit in its original packaging if it is to be stored for a long period of time, or in a safe, cool and dry place.

7. SPECIFICATIONS

VIBFIT	
Maximum power	30 W
Input	100-240V, 50/60Hz
Output	24V 1500mA
Measures	40.5x39.5x12.5 cm
Net weight / Gross weight	4.2 kg / 5.3 kg

This article is not a medical device and is not intended for the treatment or prevention of disease.

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GRUPO R. QUERALTÓ S.A.